

Operating Tables and Boards



Previous



Next



Section
Table of
Contents



Main
Table of
Contents



Product
Index



Search



WWW
Home



Contact
Us



Heated Small Animal Operating Table

- This operating table is good for general surgical work on small animals


The top is a flat, heavy gauge, biologically inert piece of stainless steel measuring 23 x 33 cm (9 x 13 in). A warming chamber located beneath this top contains a small metal sheathed 40 watt heater controlled by an energy regulator. An energy regulator controls the table top surface, varying the heating cycle when necessary. The regulator is calibrated from "full" to "off" through ten marked points. This allows adequate heating control for most operating procedures. Indicator lights are visible when the heater is on and when the regulator is switching.

Uniform heat distribution is maintained over the total surface by an insulation barrier located between the resistance mat and stainless steel top. Four stainless steel 6.5 mm (1/4 in) rods are supplied with the Table. Each 20.4 cm (8-1/16 in) rod is mounted on a simple slotted clamp with a locking screw. The slot of this clamp fits over the edge of the Table and allows the rod to be inserted anywhere along the edge of the Table or slid along the edge into position. Four tying cleats are also supplied which have a deep slot similar to the rod clamps to fit over the edge of the table and secure with a locking screw. The top of the cleat has a tapered slot that binds any tying strings securely in place. These cleats, like the rod, can be inserted at any point around the edge of the Table or can be slid along the edge into position. Folding legs under the Table allow the Table to be used in a tilted position.

Order # Product

PY2 50-1239 Heated Sm. Animal Operating Table, 115 VAC, 60 Hz

PY2 50-1247 Heated Sm. Animal Operating Table, 230 VAC, 50 Hz

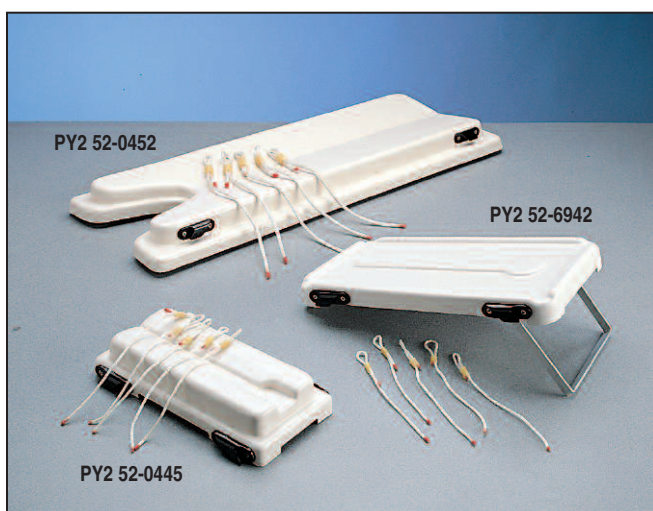

Plastic Operating Boards

- Provides excellent access to the entire length of the animal

These one-piece molded boards hold rodents and other small animals during surgery. They provide excellent access to the entire length of the animal. A quick-release system of restraint cords and locking cleats quickly and easily attach the animal to the Board. The large Rodent Board has a raised center for easily positioning the animal. Fluids drain to the sides and a folding metal leg elevates one end of the board for special positioning.

The small Rodent Board has a built-in "V" design for easy positioning of the animal. Fluids drain to the center of the Board. The Rabbit Board has a "V" design for easy positioning of the animal. Suction cups on the feet hold the unit steady. A folding leg tilts the board to a 30° angle.

Each is supplied with four cleats, one mounted at each corner, and five restraining cords for quick tie-down of animals.

Plastic Operating Boards

Order #	Board	Dimensions, L x W	Max. Weight Range	Length of Animal
PY2 59-6942	Rodent Large	40.6 x 13 cm (16 x 5-1/8 in)	Up to 500 g	35.6 cm (14 in)
PY2 52-0445	Rodent Small	27.9 x 11.4 cm (11 x 4-1/2 in)	Up to 125 g	20.3 cm (8 in)
PY2 52-0452	Rabbit	60 x 30 cm (24 x 12 in)	—	—