


ATTACHMENTS AND ACCESSORIES

2502.900.25
ORION STANDARD LARYNGEAL MIRROR STEM
WITH BULB 2.5V



7000.950.50
TIMESCO OTOSCOPE BELLOWS

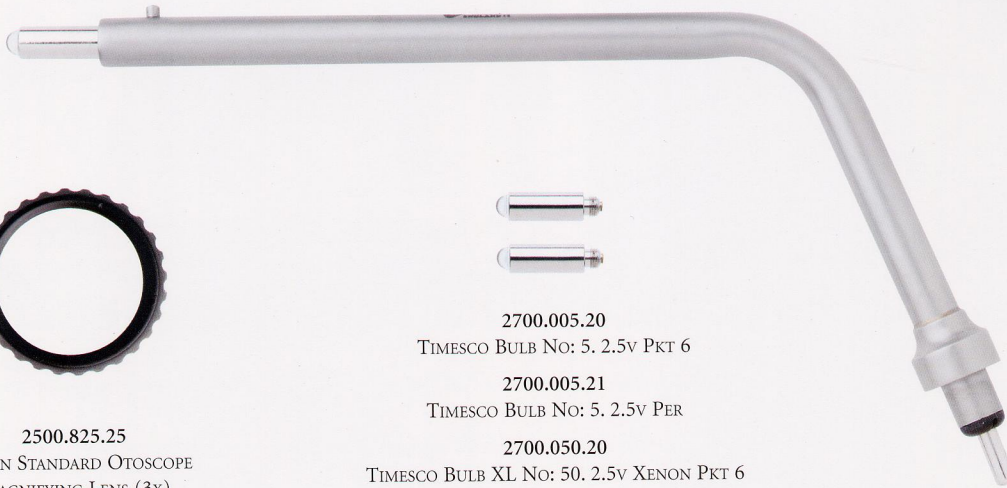

2500.825.25
ORION STANDARD OTOSCOPE
MAGNIFYING LENS (3X)



2502.926.25
ORION STANDARD LARYNGEAL STEM TONGUE
WOODEN DEPRESSOR HOLDER

2700.005.20
TIMESCO BULB NO: 5. 2.5V PKT 6

2700.005.21
TIMESCO BULB NO: 5. 2.5V PER

2700.050.20
TIMESCO BULB XL NO: 50. 2.5V XENON PKT 6

2700.050.21
TIMESCO BULB XL NO: 50. 2.5V XENON PER

2700.080.20
TIMESCO BULB XL NO: 80. 2.5V XENON PKT 6

2700.080.21
TIMESCO BULB XL NO: 80. 2.5V XENON PER

BULB NO:
5 FOR LARYNGEAL STEMS
50 FOR OTOSCOPIES
80 FOR OPHTHALMOSCOPIES




2500.920.50
ORION STANDARD OTOSCOPE NASAL SPECULUM
WITH SCREW OPENING



2502.925.25
ORION STANDARD LARYNGEAL
STEM TONGUE DEPRESSOR



2502.905.25
ORION STANDARD LARYNGEAL
STEM MIRROR NO: 3
20MM DIA

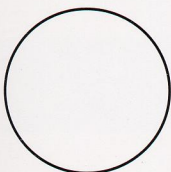

2502.905.50
ORION STANDARD LARYNGEAL
STEM MIRROR NO: 4
22MM DIA

ORION ORAL AND THROAT ATTACHMENTS

The universal fitting on the Orion handle allows various attachments that can be used to produce a greater range of Diagnostic examinations.

For Oral/Laryngeal examinations, the illuminated Orion light stem can be attached to Laryngeal mirrors (2 sizes available), a connector to use with wooden tongue

depressors or a metal Tongue depressor, whichever option is required.

Some of our components on the attachments are now being manufactured from high-density plastics, which are still durable, but have the advantage of being lightweight.

BELLOWS

For tympanic testing the Orion Otoscope has an airtight head chamber. When the Orion bellows are connected to the nipple on the otoscope head, the large bellows bulb

allows sufficient air pressure to detect movement in the tympanic membrane. The Orion bellows are made from non-cracking supple quality rubber and are free from latex.