





Introduction

The Sealsafe® PLUS range of IVC products features a PATENTED cage design, consisting of a soft silicone gasket that protrudes radially all around the cage base rim: a perfect sealing is simply achieved by placing the top on the cage body, through the radial interference between gasket and cage top.

Consequently, the vertical compression (commonly provided by the use of latches) of the silicone gasket between cage base and top is not required:

- · no latches for the best ergonomics
- no need for latches to achieve excellent containment and high protection

1/Protection for Operators and Animals

- PATENTED cage sealing design: optimal seal and protection without need for latches.
- Multiple sealed cage: silicon gasketed valves (not invasive to the cage), bottle cap, cage base and filter retainer ensure maximum sealing.
- Large Microbiological filter (VFE* ≥ 99.999987%, BFE* ≥ 99.9999937%) at cage level, ensures optimal gas exchange as well as an ideal protection also during procedures or in case of power failure.
- DOP tested Air Handling Unit's HEPA Filters provide high protection against bacteria, viruses and allergens.

*VFE: Virus Filtration Efficiency. BFE: Bacteria Filtration Efficiency.

2/Animal Welfare

- Air inlet and outlet valves in the top: no air drafts at animal level, avoiding risk of stress and heat loss, as reported in the scientific literature.
- Stand-alone Air Handling Unit: Easy, Smart and Sky Flow provide safe ventilation without transmitting vibrations to the rack, ensuring standardization of the experiment. The system is independently certified by TÜV.

3/Ergonomics

- No latches for a faster cage change, while protection is always ensured by the patented cage sealing. Ideal for Automation.
- Half length wire bar lid: easy access for a quicker cage change, minimizing the risk of musculoskeletal disorders as well as reducing the cage weight.
- Great visibility-focused design and Automatic Visual Docking Indicator with Cage Self-Locking System to help operators perform a quicker daily check.
- External self-centering water bottle: quicker bottle change without removing the cage from the rack.


4/Low Running Costs

- Ten-row racks under 2 meters high in a small footprint for maximized density: optimized floor area usage with consequent reduction of the facility running costs.
- Stand-alone air handling unit can manage up to 4 single sided racks (or 2 double-sided): energy saving with lower maintenance costs.
- ECM DC-Blowers means low power consumption and less maintenance.
- Quick daily check, easy access to the cage and quicker bottle change mean low labour costs.




Sealsafe®PLUS Mouse IVC Green Line

Protection for Operators and Animals

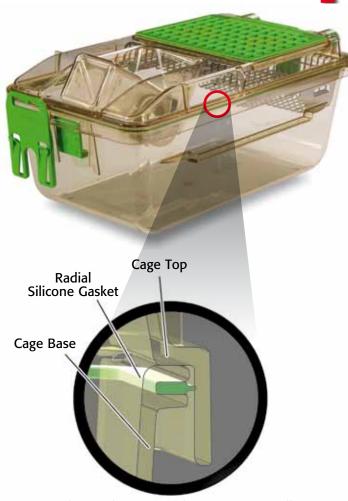

Patented cage design / Gasket-top perfect adherence: sealing with no need for vertical compression



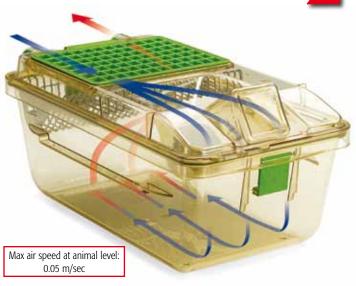

Air inlet/outlet on the top


Ergonomics








Half length lid

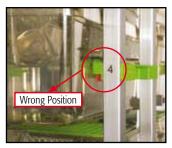

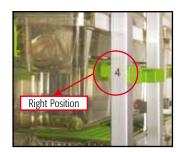

Automatic Visual Docking Indicator

Low Running Costs





- Multilinking: one AHU serves up to four racks.
- Safe for animals/operators/ environment.
- Energy-saving and Eco-compatible.
- Light and ergonomic.
- User friendly.
- High monitoring and security: Remote Alarm System (RAS) and Wireless Data Management (WDM).



TECHNICAL SHEET

GM500 cage

Floor Area:

501 cm² • 77.66 sq inch

Overall Dimensions (W x D x H): 391 x 199 x 160 (172*) mm • 15.40 x 7.83 x 6.30 (6.79*) inch *With water bottle

Plastic Type:

H-TEMP (Polysulfone) / Clear H-TEMP: 3 years warranty X-TEMP (Polyphenilsulfone): 8 years warranty

Patented Cage Body: Radial silicon gasket (occlusive seal); Simmetric (no back no front)

Half wire bar lid style

Feeder Capacity: 0.5 | • 16.91 fl oz (US) • 17.60 fl oz (UK)

Filter retainer with silicon gasket; self-centering indentation for external water bottle

Water Bottle:

260 ml • 8.80 oz (US) • 9.15 oz (UK) capacity External with o-ring for perfect sealing

Cage Filter Efficiency:

Microbiological filter

Bacterial Filtration Efficiency (BFE) ≥ 99.9999937% Virus Filtration Efficiency (VFE) ≥ 99.999987%

Sealsafe Plus rack

Cat. N.	Configuration	Specie	Overall Dimensions (L x W x H) - mm / inch
2GM54	54 cages (6 w x 9 h) Single sided	Mouse	1314 x 500 x 1850 (1949*) / 51.73 x 19.69 x 72.83 (76.73*)
2GM60	60 cages (6 w x 10 h) Single sided	Mouse	1314 x 500 x 1995 (2090*) / 51.73 x 19.69 x 78.54 (82.27*)
2GM63	63 cages (7 w x 9 h) Single sided	Mouse	1528 x 500 x 1850 (1949*) / 60.16 x 19.69 x 72.83 (76.73*)
2GM70	70 cages (7 w x 10 h) Single sided	Mouse	1528 x 500 x 1995 (2090*) / 60.16 x 19.69 x 78.54 (82.27*)
2GM72	72 cages (8 w x 9 h) Single sided	Mouse	1742 x 500 x 1850 (1949*) / 68.58 x 19.69 x 72.83 (76.73*)
2GM80	80 cages (8 w x 10 h) Single sided	Mouse	1742 x 500 x 1995 (2090*) / 68.58 x 19.69 x 78.54 (82.27*)
2GM108	108 cages (6 w x 9 h) Double sided	Mouse	1314 x 888 x 1853 (1962*) / 51.73 x 34.96 x 72.95 (77.24*)
2GM120	120 cages (6 w x 10 h) Double sided	Mouse	1314 x 888 x 1998 (2103*) / 51.73 x 34.96 x 78.67 (82.81*)
2GM126	126 cages (7 w x 9 h) Double sided	Mouse	1528 x 888 x 1853 (1962*) / 60.16 x 34.96 x 72.95 (77.24*)
2GM140	140 cages (7 w x 10 h) Double sided	Mouse	1528 x 888 x 1998 (2103*) / 60.16 x 34.96 x 78.67 (82.81*)
2GM144	144 cages (8 w x 9 h) Double sided	Mouse	1742 x 888 x 1853 (1962*) / 68.58 x 34.96 x 72.95 (77.24*)
2GM160	160 cages (8 w x 10 h) Double sided	Mouse	1742 x 888 x 1998 (2103*) / 68.58 x 34.96 x 78.67 (82.81*)

With easily removable supply plenum



Options

For Cage:

- · Cage body with grommet, for auto watering
- "Unihand" Cage body (button free in the front) for single hand procedures
- "Unihand" Cage body with grommet, for auto watering
- Full length wire bar lid with food hopper and bottle in the front, internal bottle
- Full length wire bar lid with external bottle in the front and food hopper in the back
- Flat top without bottle indentation, for auto watering
- Top with depression for disposable water pouches
- Vertical or Horizontal flip up plastic card holder
- · Raised floor grid

For Rack:

- Rack with automatic watering system
- Quick disconnecting magnetic joints (connecting AHU to rack or rack to rack)
- Flexible hoses (available in several lenghts)

Accessories

- · Mouse Balcony
- · Mouse Loft
- Mouse House

ENVIRONMENTAL SUSTAINABILITY: MINIMIZE IMPACT MAXIMIZING VALUE What Tecniplast did to become leader in sustainability:

Company:

- Quality Management System Certification according to ISO 9001 and Environmental Management System Certification according to ISO 14001.
- Certified Environmental Report according to Eni Foundation Enrico Mattei Guidelines.

Products:

- Life Cycle Assessment (LCA) according to ISO 14040 and ISO 14044.
- Carbon Footprint according to PAS 2050.
- Environmental Data Sheet according to ISO 14025.
- Contribution to credits under the LEED rating system.



















Related Products

To maximise operational efficiency and learn more about these products contact your local representative.



CS5 Evo GP changing station. Compact, ergonomic and light double-side mobile Changing Station that allows the user to easily move the cabinet inside the animal room.



Logistic trolleys offer safety and ergonomic features that enable high load-density for moving cages and bottles, crates, pallets and covers.



900 GP Rack Washer. Designed specifically to meet the throughput needs of medium to large size facilities. Maximum cleaning efficiency, safety, environmental impact and quality construction.