




2502.900.25  
ORION STANDARD LARYNGEAL MIRROR STEM  
WITH BULB 2.5V

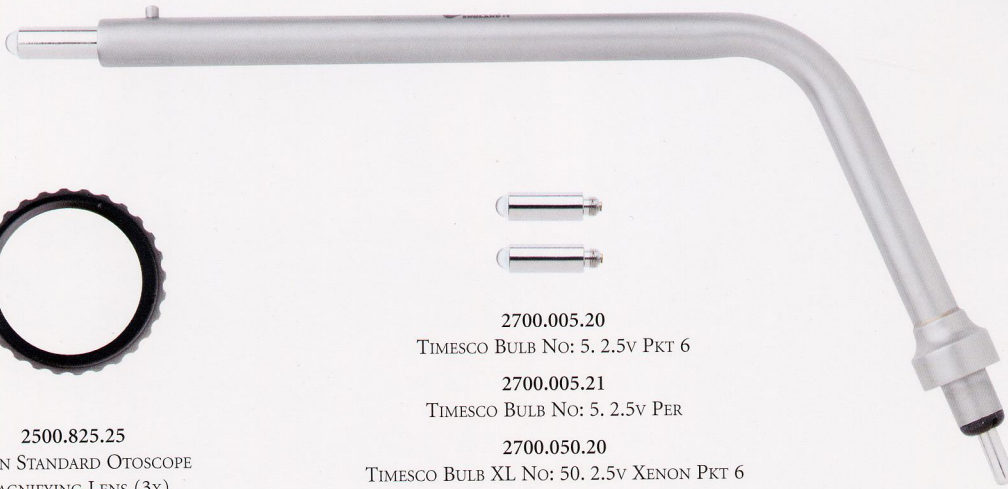

7000.950.50  
TIMESCO OTOSCOPE BELLOWS

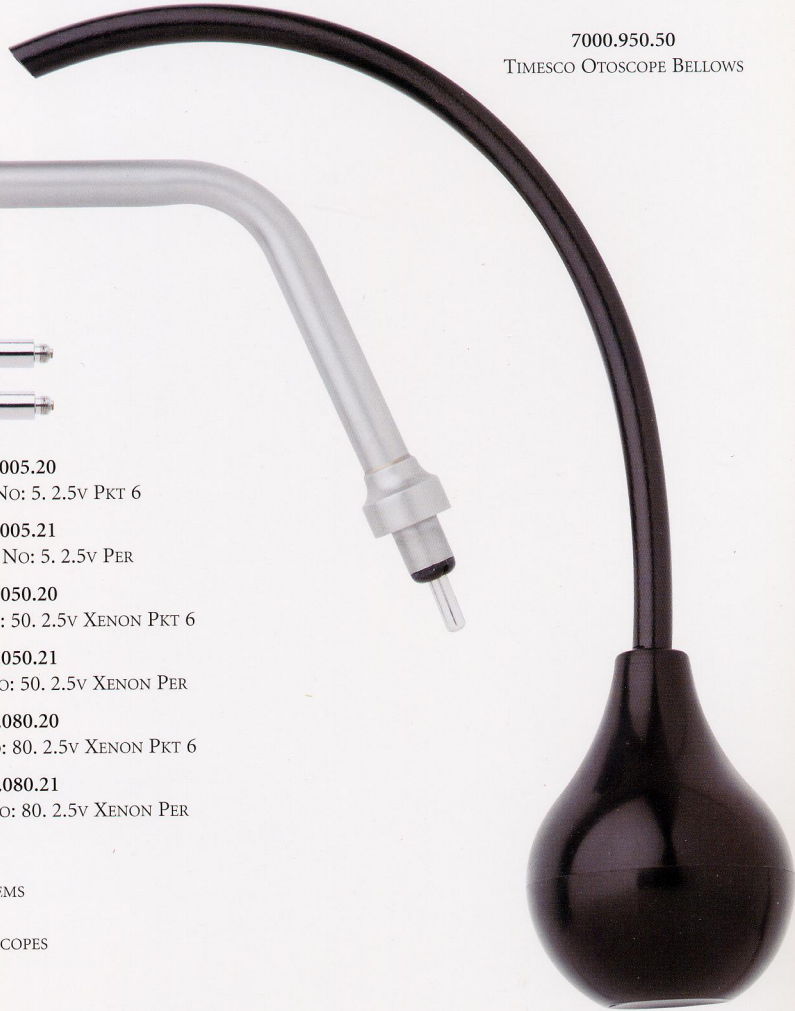

2500.825.25  
ORION STANDARD OTOSCOPE  
MAGNIFYING LENS (3X)

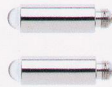

2700.005.20  
TIMESCO BULB NO: 5. 2.5V PKT 6

2700.005.21  
TIMESCO BULB NO: 5. 2.5V PER

2700.050.20  
TIMESCO BULB XL NO: 50. 2.5V XENON PKT 6

2700.050.21  
TIMESCO BULB XL NO: 50. 2.5V XENON PER

2700.080.20  
TIMESCO BULB XL NO: 80. 2.5V XENON PKT 6

2700.080.21  
TIMESCO BULB XL NO: 80. 2.5V XENON PER



2502.926.25  
ORION STANDARD LARYNGEAL STEM TONGUE  
WOODEN DEPRESSOR HOLDER

BULB NO:  
5 FOR LARYNGEAL STEMS  
50 FOR OTOSCOPIES  
80 FOR OPHTHALMOSCOPIES



2502.905.25  
ORION STANDARD LARYNGEAL  
STEM MIRROR NO: 3  
20MM DIA



2500.920.50  
ORION STANDARD OTOSCOPE NASAL SPECULUM  
WITH SCREW OPENING



2502.925.25  
ORION STANDARD LARYNGEAL  
STEM TONGUE DEPRESSOR


## ORION ORAL AND THROAT ATTACHMENTS

The universal fitting on the Orion handle allows various attachments that can be used to produce a greater range of Diagnostic examinations.

For Oral/Laryngeal examinations, the illuminated Orion light stem can be attached to Laryngeal mirrors (2 sizes available), a connector to use with wooden tongue

depressors or a metal Tongue depressor, whichever option is required.

Some of our components on the attachments are now being manufactured from high-density plastics, which are still durable, but have the advantage of being lightweight.

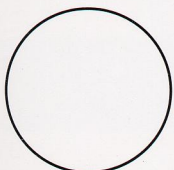

2502.905.50  
ORION STANDARD LARYNGEAL  
STEM MIRROR NO: 4  
22MM DIA

## BELLOWS

For tympanic testing the Orion Otoscope has an airtight head chamber. When the Orion bellows are connected to the nipple on the otoscope head, the large bellows bulb

allows sufficient air pressure to detect movement in the tympanic membrane. The Orion bellows are made from non-cracking supple quality rubber and are free from latex.