


RESULTS AT YOUR FINGERTIPS

- The Mini-Spirometer which provides accurate and reliable readings of PEF, FEV₁, FVC and the FEV₁: FVC ratio (as a percentage)
- Results summary and now the Volume/Time curve printed (with compatible printers)
- "BAD" blow software analysis

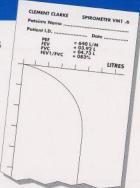

Clement Clarke

The battery powered and portable Enhanced VM1 incorporates a built-in Mini-Wright peak flow meter which enables checks to be made on the accuracy of every measurement.

Arrives complete with Normal Values guide and parallel printer connection lead.

The Enhanced VM1 is available calibrated to either ATS Scale or Wright Scale, and joins the other quality meters from Clement Clarke.

Parameters measured: PEF (L/min)

FVC (Litres)

FEV1 (Litres)

Measurement mode: Static pressure

FEV1/FVC (%)

sensing in Mini-Wright peak flow meter

Pressure Transducer: Integrated circuit

piezoresistive

Microprocessor: Z80, Clock rate 2MHz

Signal processing rate: 300Hz

Peakflow period: 10ms

Frequency Response: 100Hz @ 350L/min Sampling time: Up to 10s

Memory size: 8K RAM 4K ROM

Power Source: 1 x PP3 9v Battery Printer Compatible: Parallel, accepts 36 way IEEE 488 Centronics

Connector

Auto-switch-off time: 30 s (approx)

Range: PEF 15-800 L/min Accuracy: Better than ± 5%

(95% confidence limits)

Repeatability: Better than ± 2% Flow resistance: 12cm H₂0 at 350 L/min

Size: Including mouthpiece

210 x 120 x 70mm Weight: 420g (15oz) including

mouthpiece and battery

Flow Meter Type: Variable orifice. springloaded Piston

(Mini Wright)

Dynamic Characteristics: Critically damped at

normal profile at 350 L/Min

Catalogue

Catalogue Number

VM1 Mini-Spirometer

(inc. printer lead and Normal Values guide)

Wright Scale

3108000

ATS Scale

3108500

Number Adult disposable mouthpieces (100's) 3125030 Paediatric disposable mouthpieces (40's) 3125005 Sterilisable universal mouthpieces (15's) 3103095 Sterilisable universal mouthpieces (4's) 3103054 Pad of 50 peak flow charts 3122007



Connect With Us













